

The Beauty of Superior Anesthetic Delivery


Achieve the beauty of superior anesthetic delivery for better control of your patients' post-operative pain. The ACCUFUSER® Post-Op Pain Control Pump features the unique McKinley Saturation Catheter™ that provides thorough anesthetic coverage of the surgical site.

Designed to Saturate the Surgical Site

Available exclusively with the ACCUFUSER pump, the McKinley Saturation Catheter guides the local anesthetic through a pattern of tiny holes located around the catheter's circumference. The patented fluid conduction pattern delivers a uniform flow of medication that bathes the entire area around the wound. The result is effective post-op pain relief for your patient.

Pain Relief for Up to Five Days

The McKinley Saturation Catheter is inserted intraoperatively at the surgical site and connected to the ACCUFUSER pump flow control tubing. When the pump is activated, the catheter bathes the surgical site with continuous local anesthetic for up to five days. When the therapy is complete, the catheter removes easily, and the pump is discarded.

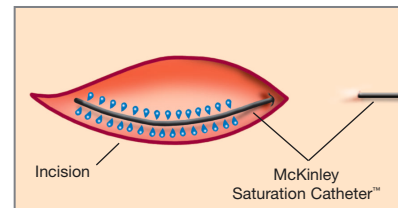


Specific to Your Surgical Needs

Designed for surgical use, the McKinley Saturation Catheter provides you these special features:

- 2 saturation lengths (2.5" & 5")
- More holes (15 & 30)
- Radiopaque
- High tensile strength

Plus, for surgeries requiring dual catheter placement, a Y-site administration set is available so you can connect two catheters to a single ACCUFUSER pump.

The beauty of personal post-op pain control™ belongs to superior anesthetic delivery using the McKinley Saturation Catheter.

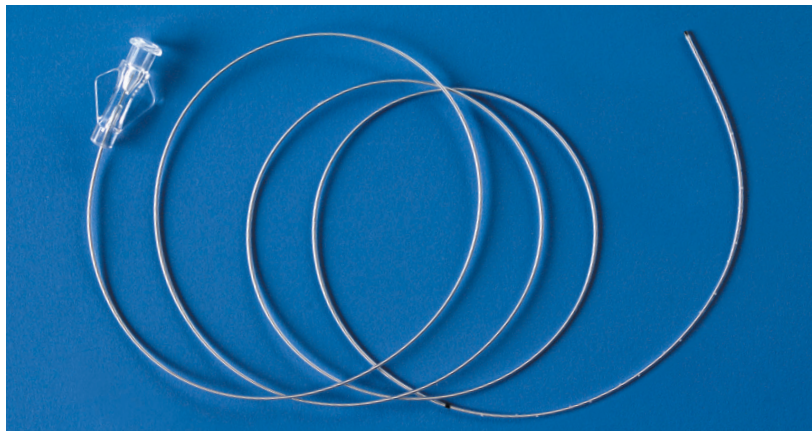

SURGICALLY-SPECIFIC FEATURES:

Patented fluid conduction pattern

Provides thorough anesthetic coverage:

- 15 holes (2.5" catheter)
- 30 holes (5" catheter)

2 saturation lengths available

- 2.5" (6.35 cm)
- 5" (12.7 cm)
- Total catheter length:
36" (91.4 cm)

Radiopaque

Verify placement and integrity

Y-site set (option)

Connect 2 catheters to 1 pump

High tensile strength

Facilitates placement

Sterile procedure tray

Contains everything you need, including:

- Catheter with fixed hub adaptor
- Introducer needle with break-away sheath
- Adhesive strip cards
- Occlusive dressing

20 gauge

MATERIALS:

Latex-free/DEHP-free

Tubing: polyether block amide

Coil: type 304 stainless steel

Fixed hub: polycarbonate



McKinley Medical

McKinley Medical, LLC

4080 Youngfield Street
Wheat Ridge, CO 80033 USA

303.420.9569

303.420.8585 (fax)

www.controlpostoppain.com

CAUTION: Federal (USA) law restricts this device to sale by or on the order of a physician.

©2005 by McKinley Medical, LLC. All rights reserved. ACCUFUSER and the McKinley Medical mountain logo are registered trademarks, and The Beauty of Personal Post-Op Pain Control and McKinley Saturation Catheter are trademarks of McKinley Medical, LLC. U.S. Patents: 6,024,724; 6,689,110B2; 6,676,643B2. PN 203663 Rev. AA 01/2005